





Base Bottom View

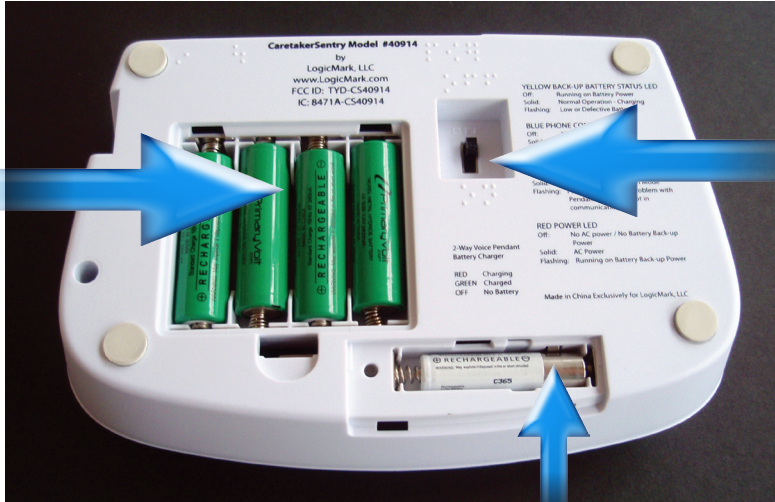

Master ON/OFF Switch

Backup Batteries
* **Note the orientation of the batteries**

Two-way Voice Pendant Battery Charging Port (for optional accessory)
* **Note the orientation of the battery**

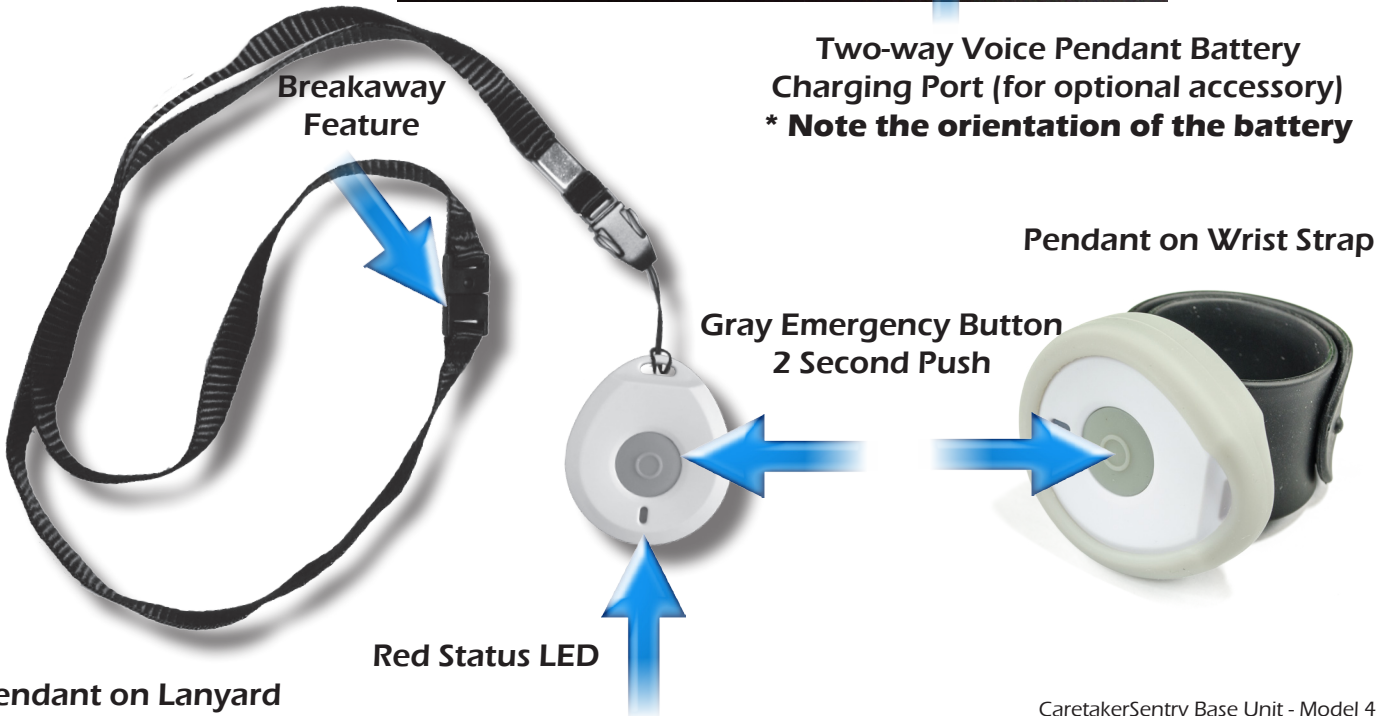

Red Status LED

Pendant on Lanyard

Pendant
Part #40915

LogicMark Disclaimers and Limited Warranty

COMMUNICATION AND RESPONSE LIMITATIONS: Purchaser acknowledges that signals which are transmitted over telephone lines, or other modes of communication pass through communication networks wholly beyond the control of LogicMark LLC and are not maintained by LogicMark LLC, and, therefore, LogicMark LLC shall not be responsible for any equipment or communication failure which prevents transmission signals from reaching your contact list including emergency 911 operators or damages arising therefrom. **Purchaser acknowledges that LogicMark LLC provides no response to the System's equipment.** The equipment is designed to communicate with the a central station of your choice and LogicMark LLC is not and shall not be responsible for ambulance, police or other emergency response time or that any response will be provided by the central monitoring station.

TESTING AND SERVICE OF THIS EQUIPMENT: The equipment, once installed, are in the exclusive possession and control of the Purchaser, and it is Purchaser's sole responsibility to test the operation of equipment and request warranty service if the equipment is under warranty.

PURCHASER'S EXCLUSIVE REMEDY: Purchaser's exclusive remedy for LogicMark LLC's default hereunder is to require LogicMark LLC to repair or replace, at LogicMark LLC's option, any equipment or part of the personal emergency alert system which is non-operational during LogicMark LLC's warranty period.

LIMITATION OF LIABILITY: This equipment is not designed or guaranteed to prevent any loss or injury. This Limited Warranty and Disclaimer of Liability constitutes the terms of sale and use of the equipment, and if there should arise any liability on the part of LogicMark LLC as a result of any cause whatsoever, regardless of whether or not such loss, damage, or personal injury was caused by or contributed to by LogicMark LLC's negligence to any degree or failure to perform any obligation or strict products liability, such liability will be limited to an amount paid by the Purchaser to LogicMark LLC for the product, or to the sum of \$350.00, whichever is greater.

For warranty information, see your dealer for details.

CARETAKERSENTRY


Quick Install Guide

Instructions for Installing Your New Personal Emergency Response System

Base Top View
Part #40914



Base Back View

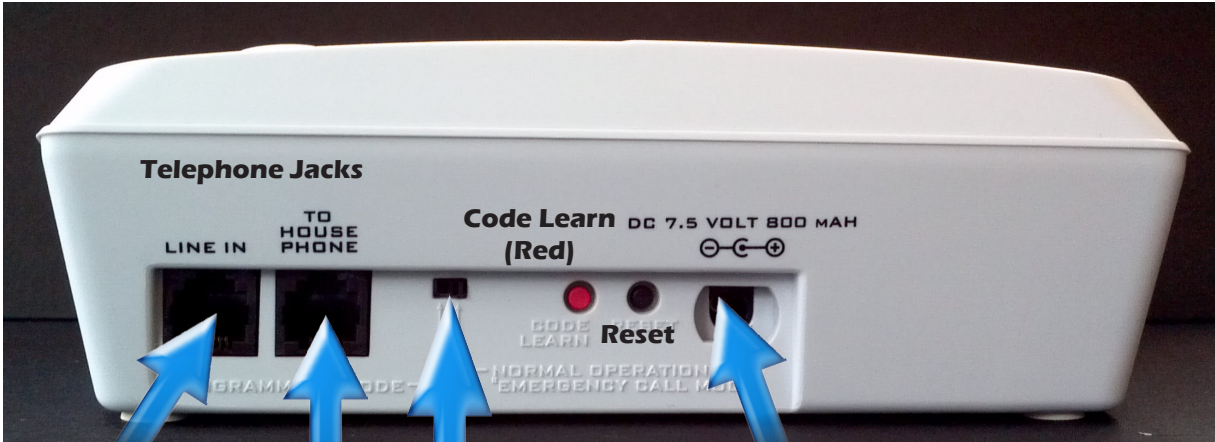

Line IN

To House Phone

Power DC 7.5
800 mAh

Programming Switch
Programming Mode "Emergency Call Mode"
Normal Operation

Base Side View



Cellular Accessory Jack

CARETAKERSENTRY Step-by-Step Installation

1

Confirm Accessories in Kit

4 Rechargeable AA Batteries



Telephone Cord



AC Adapter



Lanyard




Wrist Strap
(Shown with basic pendant)



Model #40911

3

Plug Telephone Line and Power Cord into Wall Jacks

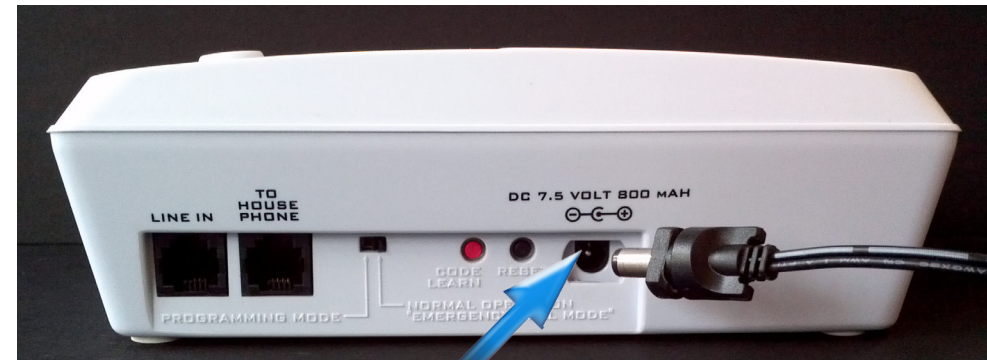


Plug the AC Adapter power cord into the power jack and twist 1/4 turn clockwise to lock.



2

Place Batteries in Unit

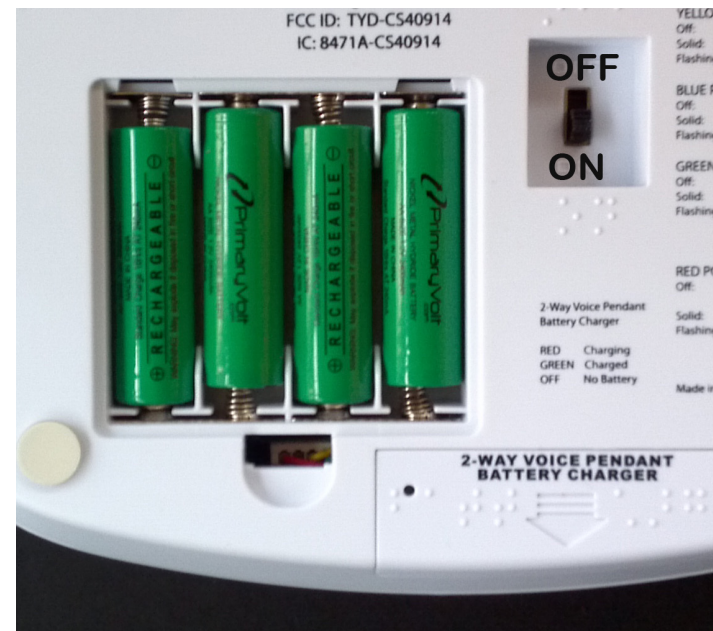

Place Green Rechargeable Batteries in the compartment on the bottom of the base unit. **Make sure to turn ON/OFF switch to ON. You will hear the announcement "Running on Battery Power." Unit will do so until the AC Adapter power supply is plugged in.**

4

Test the Unit

Performing a System Test

- PUSH and HOLD the CANCEL button for **8-10 seconds** to enter **"System Check"**
- Push GRAY EMERGENCY button for **2 seconds** on Pendant (**RED LED on the pendant should turn ON for 1 second**)
- Listen for voice prompt from base unit:
 - "All Systems are OK – Battery is OK" is heard if the system is correctly installed.
 - If any other voice prompt is heard, the unit will describe the problem
- Push and release the CANCEL button for 1 second to exit Test Mode (unit will automatically exit in 3 minutes)
- We recommend that you perform this System Test weekly**

Place call to Emergency Operator:

- Push GRAY EMERGENCY button for 2 seconds on Pendant (RED LED on the pendant should turn ON for 1 second)
- You should hear a dial tone, dialing and then communication tones between the CaretakerSentry and central station. Several seconds later the operator should answer the call
- Talk normally to the base unit – just like speaker phone
- At the end of the conversation, the CaretakerAlert should automatically hang-up – if not, push the CANCEL button for 1 second

